

EQUI RECOVER

Instruction Manual Of Cold Therapy System

[Name]: Cold Therapy System

[Description]: The product consists of cooler, tube and polyurethane coated nylon.

[Intended Use]: The Cold Therapy System is intended to be used for physical antipyretic and cold therapy on specific parts of human body surface.

[Use Instruction]:

- 1. Prepare cooler:** connect the blue tube to the cooler. Add water to line inside the cooler, then add ice. For 6 to 8 hour treatment, add ice to top of cooler. Lay the insulation disk on top of ice. Attach cooler lid snugly. Allow 5 minutes with occasional shaking to chill water.
- 2. Apply the empty cuff:** Use a medical isolation pad between your desired area and the cuff, Position cuff on thigh with quick-disconnect pointing down (Fig.1). Thread straps through white rings and secure snugly (Fig.1). Cuff can be positioned on front or back of thigh.
- 3. Fill cuff:** Connect blue tube to cuff. Open air vent on cooler lid and raise cooler no more than 15.7 inches (40cm) above the cuff (Fig.2) (if the cooler is raised more than 15.7 inches (40cm), the cuff will increase in volume and weight.) Hold raise cooler for about 30 seconds while cuff fills. Close air vent on cooler lid. the cooler can now be disconnected from cuff by pressing the tab on the quick-disconnect while cooler is raised.
- 4. Recharge water:** At first, recharge the water in the cuff after 15 or 30 minutes, then once an hour or as needed. Reconnect the blue tube to cuff. lower the cooler below cuff, and the warmed water will drain from the cuff into the cooler. Allow a minute or two for the water to mix with the ice and recharge, then raise the cooler and repeat the filling process (see step 3).

[Warning]: 1. This product does not directly contact the surface of horse or human skin and can be used repeatedly.


2. Wipe the surface of the cuff with alcohol every week for disinfection. Connect tube and cooler to flush with water.
3. After disinfection or cleaning and drying, store them flat for later use.
4. After use, completely drain water from cuff, tube, and cooler. (To drain tube, elevate tube while pressing tip of quick-disconnect.) Store cooler with top off to allow drying.

[Contraindication]: Use with caution if allergic to the product.

[Storage]: This product should be stored in a temperature range is -10 °C to 30 °C, relative humidity is not more than 80%, well ventilated room.

[Service Life]: Three years

class: Class I

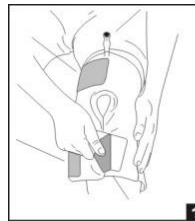

EC REP

MedNet EC-REP GmbH
Borkstrasse 10,
48163 Muenster, Germany

Equi recover
Opperuit 360
2941 AR Lekkerkerk
Netherlands

size: 141mm*210mm